


LymphFlow

GO 



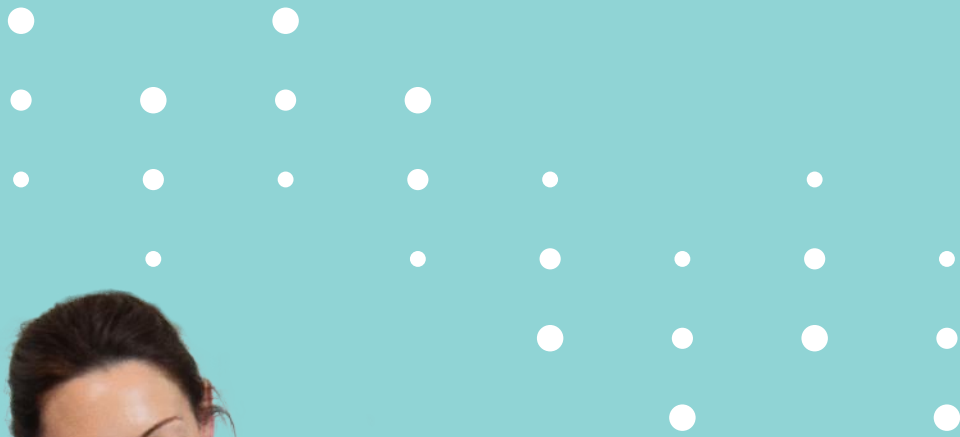

Go with the flow

The Haddenham LymphFlow Go is an all-in-one pneumatic compression device designed to aid in reducing swelling related to primary or secondary lymphoedema.

Featuring a cord free design, this innovative device combines the power of a standalone pump with a compact, portable form.



Using the latest evidence


The Lymphflow Go mimics Manual Lymphatic Drainage (MLD) techniques to stimulate the lymphatic system and direct fluid to functional lymph nodes.

12 Sequentially Inflatable Chambers

Treatment modes designed by lymphoedema specialists to promote lymphatic drainage and relieve symptoms.

Therapeutic Pressures

Tailor your treatment with a range of pressure settings to optimise results - ranging from 40 - 80mmHg.

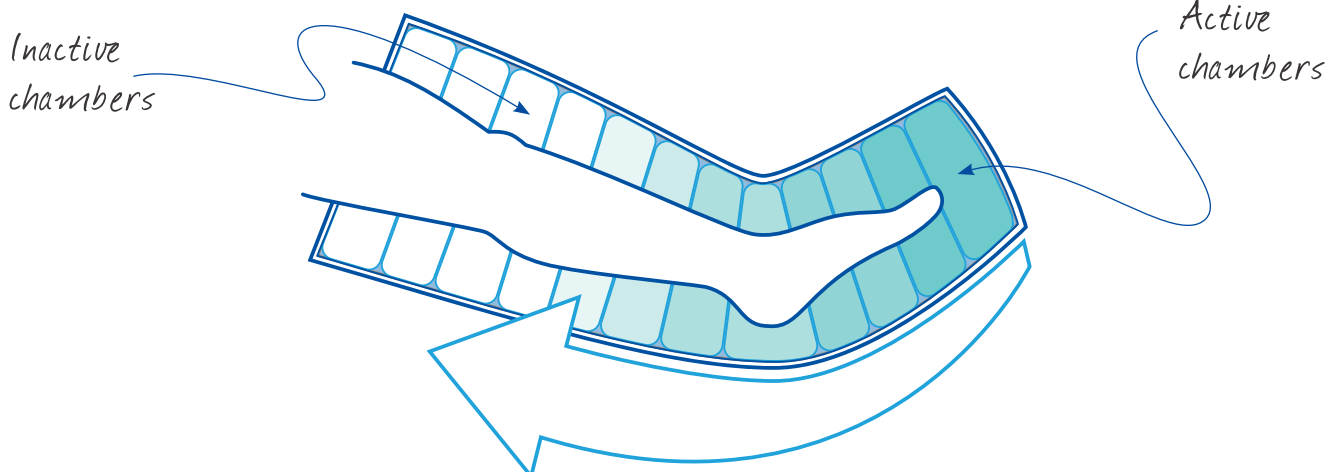

Save time for quality time.

The Lymphflow Go can seamlessly integrate into your daily routine, making it easy to incorporate therapy sessions into your everyday life. Perfect for storytime or just some me time!

Rechargeable

Up to 5 treatment cycles on a single charge.

Efficient

Treatment times ranging from only 20-45 minutes*

Convenient

Setup anywhere you can elevate your leg.



*Treatment time can vary depending on leg size, treatment mode, and pressure level.

Compact & cord free

The unique all-in-one design includes an integral pump module, which means convenient cordless usage and a lightweight treatment solution on the go.



Compact

The device packs neatly into its own carry bag, which fits within a car boot or a suitcase*

Cordless

Cordless operation means a simple setup for treatment anywhere, at anytime.

* When flying Lymphflow Go cannot be included in checked luggage due to lithium battery.

Supports self-management

“As a therapist, I’m keen to empower my patients and help them manage their lymphoedema through the challenges of everyday life.” Karen Cooper, a Manual Lymphatic Drainage practitioner with over 20 years of experience, reports that the Lymphflow Go has provided her patients with a practical self-care plan, enabling them to take control of their lymphoedema symptoms.

“LymphFlow go is an easy, affordable way to maintain the good results we achieve in clinic.”

Karen Cooper

Manual Lymphatic Drainage practitioner

Simple sizing

The 12 air chambers in the Lymphflow Go adapt to accomodate a broad range of shapes and sizes.

Take a circumference measurement at the top of your affected thigh to see which will be the best fit for you - Regular or Large.



Controls at your fingertips

Everything you need to control the Lymphflow Go is within easy reach - from pausing treatment to viewing battery level.



Mode

Set the treatment mode
Lymphatic or Sequential.

See reverse for more information



Start/Pause

Pause treatment
mid-cycle and resume
where you left off.



Level

Toggle the pressure:

Low (40mmHg)

Medium (60mmHg)

High (80mmHg)



Treatment modes



The Lymphflow Go has 2 modes which can be used in the treatment of primary or secondary lymphoedema.



Lymphatic Mode

for use where the lymphatic drainage pathway has been disrupted, such as following lymph node removal or radiotherapy to the proximal/regional nodes.



Sequential Mode

is designed for use where the proximal/regional lymph node drainage is intact, including primary lymphoedema and secondary lymphoedema without nodal disruption.

The Lymphflow Go is intended for use by medical professionals and adult patients who are under medical supervision. If unsure please consult with your therapist or doctor before purchasing.

Not suitable if you suffer from

Severe heart problems, Recent serious vein problems, Serious blood flow issues in the legs, Infections or inflammation in your skin or limbs, Situations where increased fluid return is not safe.

Ordering

The Lymphflow Go is available exclusively via Haddenham Healthcare's patient direct website Lymphshop www.lymphshop.com



LYMPHSHOP

Easy referral

Tell us your therapist's details and we'll handle the referral process.

Spread the cost

Pay for your Lymphflow Go in installments.



Patient-friendly
Compression

Call 01844 208842
Email sales@hadhealth.com
www.haddenham.com